



AMSCO® Flexmatic® Scrub Sink

Convenient. Durable. Clean.

STERIS Heavy Gauge Stainless Steel Scrub Sinks


## Designed for Safety

- STERIS' expansive design reduces undesired physical contact with the sink and minimizes risk of contamination
  - Minimal welds help prevent bacteria from collecting
- **Integrated eye-wash station** - faucet attachment provides convenient OR access with an emergency eye wash station
  - When water is turned on with knee panel, one pull of the center knob releases both orange safety caps to provide instant eye flushing\*
  - When knob is not activated, faucet functions as a scrub sink
- **Optimum water temperature** - the built-in thermostatically controlled mixing valve provides water temperatures of up to 115° F (46° C) to prevent scalding



Optional laminar flow head with optional eye wash.

Choice of rose spray or laminar flow

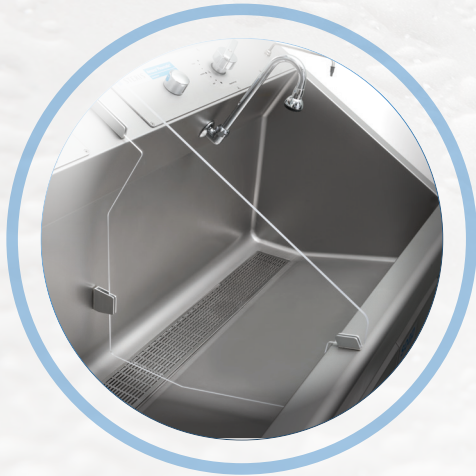

Optional divider prevents cross contamination



\*Faucet-mounted eye wash stations (regardless of manufacturer) require two motions to operate them (turn water on, pull knob) and thus do not meet the provisions of ANSI Z358.1-2004; these are intended as supplemental installations to other dedicated eye washing stations.

AMSCO® Flexmatic® Scrub Sink

## Built for Your OR

- Easy cleaning – continuous welded compartments and covered compartments constructed entirely of stainless steel
- Timed scrubbing - optional digital timer ensures adequate scrubbing time
  - Optional infrared sensor turns water on and off
  - Upgradeable knee-operated internal soap dispenser
  - **Available models:** 1, 2 or 3 station



Digital timer



Optional infrared sensor

### Convenience with the push of a knee.

Efficient hand washing starts with hands free water and soap control. Knee-activated panel turns the water on and off and dispenses the right amount of soap every time.



AMSCO® Flexmatic®  
Double Scrub Sink

# Technical Specifications

Models	AMSCO Flexmatic Scrub Sink		
	1-station	2-station	3-station
Dimensions (W x D x H)	32 x 27 9/32 x 51 1/2" (813 x 692 x 1308 mm)	64 x 27 9/32 x 51 1/2" (1626 x 692 x 1308 mm)	96 x 27 9/32 x 51 1/2" (2438 x 692 x 1308 mm)
Laminar flow faucet head	optional	optional; both faucets would have same style head	optional; all faucets would have same style head
Eye wash station attachment	optional	optional; can be attached to either left or right faucet	optional; can be attached to any of the three faucets
Infrared water control	optional	optional	optional
Foot pump soap dispenser	optional	optional; recommend one for each station	optional; recommend one for each station
Knee-operated internal soap dispenser	optional	optional	optional
Digital timer	optional	optional	optional
Divider screen	N/A	optional	optional
Temperature safety	≤ 115°F (46°C)	≤ 115°F (46°C)	≤ 115°F (46°C)

For more information, contact your STERIS representative or visit [www.steris.com](http://www.steris.com).



STERIS Corporation  
5960 Heisley Road  
Mentor, OH 44060-1834 · USA  
800-548-4873  
[www.steris.com](http://www.steris.com)